

AmbiFi allows media resource to be attached to items, this helps in providing additional rich content to support the narrative content provided in the labels which Mira will speak.

The activity has a special media item called the Title Image, this is the image shown on the card of the activity when it shows up in a users inventory.

The screenshot displays the AmbiFi Training interface. On the left is a sidebar with a top bar labeled 'AmbiFi Training' and a menu containing 'List 1', 'Item Types', 'Media Resources', and 'Media Item'. The main area on the right is a form for creating or editing an activity. The form has several sections: 'Miscellaneous' with a 'Description' field; 'Drive Root Path' with a text field containing '/private/activities/ ambifi-training'; 'Title Image' with a button labeled 'Image' (highlighted by a red arrow); 'Media' with dropdown menus for 'Background', 'Image', 'Video', 'Audio', and 'PDF'; 'General' with a 'Tags' field and 'Creator' set to 'Admin User'; and 'Basic Options' at the bottom.

Click on the "Image" button down arrow to select where to get the image from

There are three options:

- Upload a new one (default when clicking the Image button)
- External URL - if you have a URL to an image
- From Drive - if you previously uploaded an image and want to use it again

The screenshot displays the Ambifi Training interface. On the left is a sidebar with a top bar labeled "Ambifi Training" and a list of items: "List 1", "Item Types", "Media Resources", and "Media Item". The main area on the right contains several sections: "Miscellaneous" with a "Description" field; "Drive Root Path" with a text input showing "/private/activities/ ambifi-training"; "Title Image" with a dropdown menu currently set to "Image". This dropdown menu is open, showing four options: "Upload", "From External URL", "From Drive", and "Paste Resource". A red arrow points to the "Upload" option. Below the dropdown are buttons for "Background", "Image", "Video", "Audio", and "PDF". Further down is the "General" section with "Tags" and "Creator" fields, and at the bottom is the "Basic Options" section.

Once the image is selected you will see a status spinner and text (Uploading image...), once the image is uploaded you will see its preview.

The screenshot displays the AmbiFi Training web application interface. On the left is a sidebar with a blue header containing a question mark icon and the text "AmbiFi Training". Below the header are four menu items: "List 1" (with a minus icon), "Item Types" (with a plus icon), "Media Resources" (with a minus icon), and "Media Item" (with a checkmark icon). The main content area on the right is titled "Title Image" and features an "Image" button. Below this, there is a section labeled "Image" containing a "From AmbiFi Drive" button and a red "X" icon. A file path is displayed: "private/activities/ambi-fi-training/photo-of-woman-writ", followed by a "Browse" button and a dropdown arrow. A large red arrow points to the image preview, which shows a woman in a white shirt writing on a chalkboard while two men observe. At the bottom of the preview is a "View In New Tab" button.

Items also have a media toolbar, this toolbar allows to add certain media types to an item, lets begin with an image, click the Image button.

Note: On mobile devices when clicking the media buttons for image, video or audio you will be presented with options to either capture/record media or select existing media from your device.

The screenshot displays the Ambifi Training interface. On the left is a sidebar with a top toolbar containing icons for various actions. The sidebar menu includes 'Ambifi Training', 'List 1', 'Item Types', 'Media Resources', and 'Media Item With Image' (which is highlighted with a blue bar and a checkmark). The main content area on the right is titled 'Media Item With Image' and contains several sections: a 'Preview' section showing 'Media Item With Image', a 'Label 2' section with an input field containing '(Optional: Label 2)' and a 'Don't Speak Label 2' checkbox, a 'Media' section with a toolbar containing 'Background', 'Image', and 'Audio' buttons (with a red arrow pointing to the 'Image' button), and a 'PDF' button. Below these are expandable sections for 'Comments', 'Basic Options', 'Styling', 'Print Options', and 'Advanced Options' (which is highlighted in green).

As previously shown for the title image, once the image is uploaded it is presented as a preview.

The screenshot displays the 'Media Item With Image' configuration interface. On the left is a sidebar with a top toolbar and a list of items: 'Ambifi Training', 'List 1', 'Item Types', 'Media Resources', and 'Media Item With Image' (which is selected and highlighted in blue). The main area contains the following sections:

- Media Item With Image**: A text input field.
- Don't Speak Label 1**: A checkbox.
- Preview**: A section showing the 'Media Item With Image'.
- Label 2**: A section with an '(Optional: Label 2)' text input field and a 'Don't Speak Label 2' checkbox.
- Media**: A section with a toolbar for selecting media types: Background, Image, Video, Audio, and PDF.
- Image**: A section with a 'From Ambifi Drive' button and a file path 'private/activities/ambifi-training/hand-hygiene.jpg'. A 'Browse' button is also present.
- Preview**: A large preview of the 'How To Wash Your Hands' infographic, which is a 2x5 grid of numbered steps (01-10) illustrating proper handwashing technique. A red arrow points to this preview image.

To upload a video repeat the previous steps but click the "Video" button instead.

Note: When clicking the video or other media upload buttons on a PC the file browser will be limited to files that are supported. You can't upload a PDF when clicking the "Video" button.

The screenshot displays the Ambifi Training interface. On the left, a sidebar lists several options: 'AmbiFi Training', 'List 1', 'Item Types', 'Media Resources', 'Media Item With Image', and 'Media Item With Video'. The 'Media Item With Video' option is highlighted with a blue background and a white checkmark. The main content area on the right is titled 'Media Item With Video' and contains several sections. The 'Label 1' section has a checkbox for 'Don't Speak Label 1'. The 'Label 2' section has a text input field with the placeholder '(Optional: Label 2)' and a checkbox for 'Don't Speak Label 2'. The 'Media' section features a row of buttons: 'Background', 'Image', 'Video', and 'PDF'. A red arrow points to the 'Video' button. Below this, the 'Video' section shows a dropdown menu with 'From AmbiFi Drive' selected. A red 'X' button is to the right. Below the dropdown, a file path is shown: 'private/activities/ambifi-training/ambifi-training', with a 'Browse' button to its right. A red arrow points to the 'Browse' button. Below the file path, a video player is shown with the title 'AmbiFi - Getting Started' and a progress bar indicating 0:03 / 1:48. A 'View In New Tab' button is at the bottom.

The same steps again to upload a PDF

The screenshot shows the Ambifi Training interface. On the left sidebar, under 'Media Resources', there are three options: 'Media Item With Image', 'Media Item With Video', and 'Media Item With PDF'. The 'Media Item With PDF' option is highlighted with a blue background and a checkmark. A red arrow points from this option to the 'Media' section on the right. In the 'Media' section, there are buttons for 'Background', 'Image', 'Video', and 'PDF'. A red arrow points from the 'PDF' button to the 'PDF' section below. The 'PDF' section has a dropdown menu set to 'From Ambifi Drive' and a 'Browse' button. Below this, a preview of a PDF document is shown. The document is the 'REPROCESSING MANUAL' for the 'OLYMPUS EVIS EXERA II' GASTROINTESTINAL VIDEOSCOPE. It lists various models and accessories.

(Optional: Label 2)

☐ Don't Speak Label 2

Media

Background Image Video PDF

PDF

From Ambifi Drive

private/activities/ambifi-training/Olympus-GI-Endoscc Browse

OLYMPUS REPROCESSING MANUAL

INSTRUCTIONS

EVIS EXERA II

EVIS EXERA II GASTROINTESTINAL VIDEOSCOPE

OLYMPUS GIF TYPE N180
OLYMPUS GIF TYPE XP180N
OLYMPUS GIF TYPE Q180
OLYMPUS GIF TYPE H180
OLYMPUS GIF TYPE H180J

EVIS EXERA II COLONOVIDEOSCOPE


OLYMPUS CF TYPE Q180AL/I
OLYMPUS CF TYPE H180AL/I
OLYMPUS CF TYPE H180DL/I
OLYMPUS PCF TYPE Q180AL/I
OLYMPUS PCF TYPE H180AL/I

Accessories:

- Water resistant cap (MH-553, MAJ-942)
- Biopsy valve (MB-338)
- Air/water valve (MH-438)
- Auxiliary water tube (MAJ-835)
- Suction cleaning adapter (MH-456)
- Channel cleaning brush (BW-307)
- Chain for water-resistant cap (MAJ-1519)
- Suction valve (MH-443)
- Mouthpiece (MB-142, MA-474)
- Auxiliary water inlet cap (MAJ-215)
- Channel plug (MH-546)
- Injection tube (MH-946)

All resources that you upload end up in the AmbiFi drive, think of this as being something like the Google Drive, Microsoft One Drive or similar.

All files will be stored securely there and can now be organized, renamed and updated in a central place. This also allows others to reuse all the media you upload.

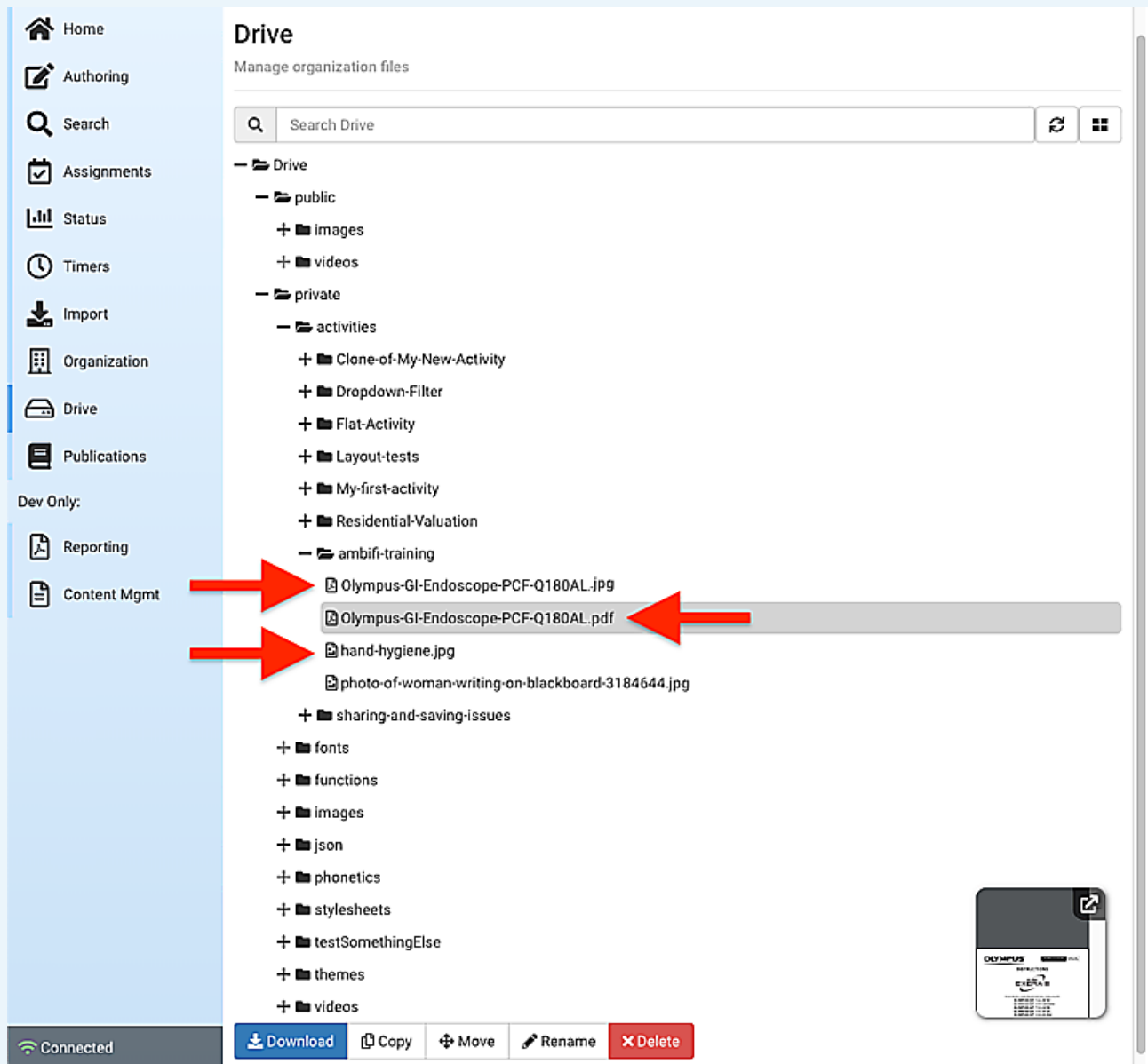


Lets now create a final item with all three media assets we already uploaded (image, video and PDF).

Select the down arrow next to image and select "From Drive" (remember we already uploaded the media)

The screenshot displays the Ambifi Training interface. On the left, a sidebar lists several items: 'AmbiFi Training', 'List 1', 'Item Types', 'Media Resources', 'Media Item With Image', 'Media Item With Video', 'Media Item With PDF', and 'Media Item With All Previous Resources' (highlighted in blue). The main area on the right is titled 'Media Item With All Previous Resources'. It includes a 'Preview' section showing the item name. Below this, there are sections for 'Label 2' (with an optional label input) and 'Media'. The 'Media' section features a row of media type selectors: 'Background', 'Image', and 'Audio'. The 'Image' selector is open, showing a dropdown menu with options: 'Upload', 'From External URL', 'From Drive' (highlighted with a red arrow), and 'Paste Resource' (also highlighted with a red arrow). Below the media section are expandable sections for 'Comments', 'Basic Options', 'Styling', 'Print Options', and 'Advanced Options'.

Once you selected the "From Drive" option make sure to "Browse" the drive and select the image.

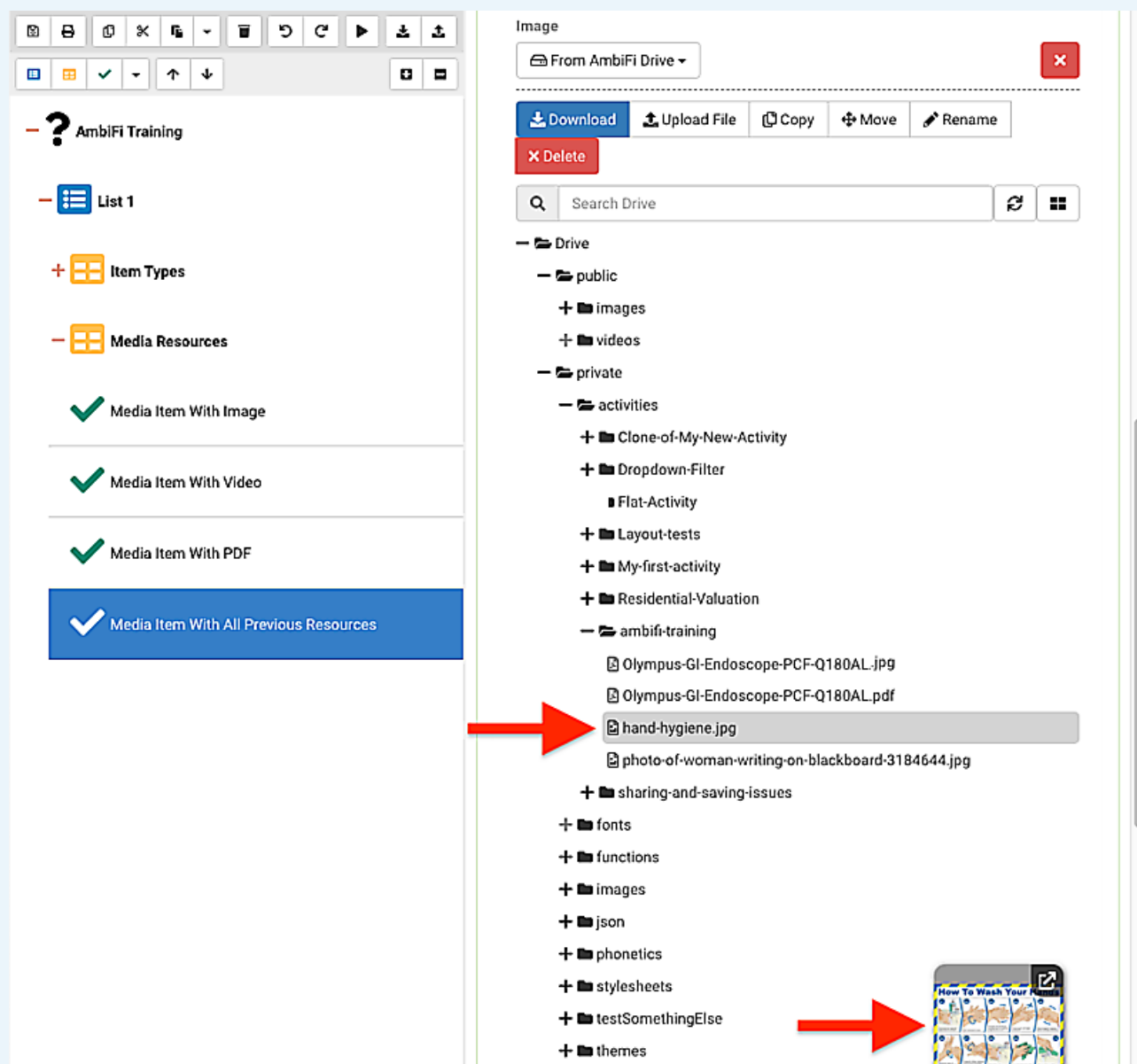

Repeat the same steps for video and PDF to embed all three in a single item.

The screenshot displays the AmbiFi Training interface. On the left is a sidebar with a menu:

- AmbiFi Training
- List 1
- Item Types
- Media Resources
 - Media Item With Image
 - Media Item With Video
 - Media Item With PDF
 - Media Item With All Previous Resources

Red arrows point from the 'Media Item With Video' and 'Media Item With All Previous Resources' options to the main content area. The main content area is divided into two sections:

Video Section:

- URL: `private/activities/ambi-fi-training/ambi-fi-training-`
- Video Player: Shows a video titled "AmbiFi - Getting Started" with a progress bar at 0:03 / 1:48.
- Button: "View In New Tab"

PDF Section:

- PDF Source: "From AmbiFi Drive"
- URL: `private/activities/ambi-fi-training/Olympus-GI-Endoscc`
- PDF Document: A document titled "OLYMPUS REPROCESSING MANUAL" with the following content:
INSTRUCTIONS
EVIS EXERA II
EVIS EXERA II GASTROINTESTINAL VIDEOSCOPE
OLYMPUS GIF TYPE N180
OLYMPUS GIF TYPE XP180N
OLYMPUS GIF TYPE Q180
OLYMPUS GIF TYPE H180
OLYMPUS GIF TYPE H180J
EVIS EXERA II COLONOVIDEOSCOPE
OLYMPUS CF TYPE Q180AL/I
OLYMPUS CF TYPE H180AL/I
OLYMPUS CF TYPE H180DL/I
OLYMPUS PCF TYPE Q180AL/I
OLYMPUS PCF TYPE H180AL/I

And now preview the activity to see all media in a single item

